

### 29.6.3 Chemistry Paper 3 (233/3)

1 You are provided with:

- acid A labelled solution A;
- 2.0 M sodium hydroxide solution labelled solution B;
- Solution C containing 25.0 g per litre of an alkanolic acid.

You are required to:

- (a) prepare a dilute solution of sodium hydroxide, solution B.
- (b) determine the:
  - (i) molar mass of the alkanolic acid
  - (ii) reaction ratio between sodium hydroxide and acid A.

#### Procedure 1

Using a pipette and a **pipette filler**, place 25.0 cm<sup>3</sup> of solution B into a 250.0 ml volumetric flask. Add about 200 cm<sup>3</sup> of distilled water. Shake well. Add more distilled water to make up to the mark. Label this solution D. Retain the remaining solution B for use in procedure II.

Fill a burette with solution C. Using a clean pipette and a **pipette filler**, place 25.0 cm<sup>3</sup> of solution D into a 250 ml conical flask. Add two drops of phenolphthalein indicator and titrate with solution C. Record your results in **table 1**. Repeat the titration two more times and complete the table.

Table 1	I	II	III
Final burette reading			
Initial burette reading			
Volume of solution C (cm <sup>3</sup> ) added			

(4 marks)

Determine the:

- (i) average volume of solution C used; (1 mark)
- (ii) concentration of solution D in moles per litre; (1 mark)
- (iii) concentration of the alkanolic acid in solution C in moles per litre (1 mole of the acid reacts with 3 moles of the base); (1 mark)
- (iv) molar mass of the alkanolic acid. (1 mark)

#### Procedure II


Fill a **clean** burette with solution A. Place 5 cm<sup>3</sup> of solution A into a 100 ml beaker. Measure the initial temperature of solution A in the beaker and record it in **table II**. Using a 10 ml or a 100 ml measuring cylinder, measure 25 cm<sup>3</sup> of solution B. Add it to solution A in the beaker and immediately stir the mixture with the thermometer. Record the maximum temperature reached in **table II**. Repeat the experiment with other sets of volumes of solutions A and B and complete the table.

**Table II**

Volume of solution A (cm <sup>3</sup> )	5	9	13	17	21	25
Volume of solution B (cm <sup>3</sup> )	25	21	17	13	9	5
Maximum temperature (°C)						
Initial temperature (°C)						
Change in temperature, $\Delta T$						

(6 marks)

- (a) On the grid provided; plot a graph of  $\Delta T$  (Vertical axis) against the volume of solution A. (3 marks)

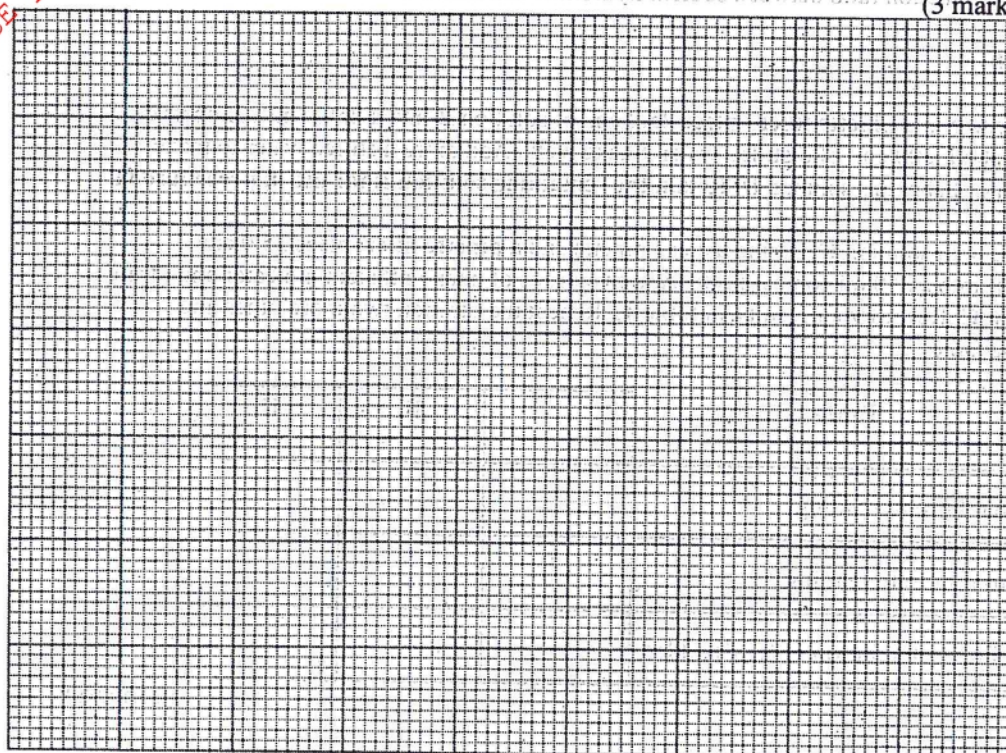

- (b) From the graph, determine the volume of solution A which gave the maximum change in temperature. (1 mark)
- (c) Determine the volume of solution B that reacted with the volume of solution A in (b) above. (1 mark)
- (d) Calculate the:
- ratio between the volumes of solutions A and B that neutralised one another; (1 mark)
  - concentration in moles per litre of the acid in solution A. (Assume that the volume ratio is the same as the mole ratio). (1 mark)



2

You are provided with solids E, F and G.

Carry out the tests below and write your observations and inferences in the spaces provided.

- (a) Place all of solid E in a boiling tube. Add 20 cm<sup>3</sup> of distilled water and shake until all the solid dissolves. Label this as solution E.

- (i) To about 2 cm<sup>3</sup> of solution E in a test-tube, add 4 drops of 2 M sulphuric (VI) acid.

OBSERVATIONS

INFERENCES

(1 mark)

(2 marks)

- (ii) To about 2 cm<sup>3</sup> of solution E in a test-tube, add 2 M sodium hydroxide dropwise until in excess.

OBSERVATIONS

INFERENCES

(1 mark)

(1 mark)

- (iii) Place one half of solid F in a test-tube. Add 2 cm<sup>3</sup> of distilled water and shake well. Add 4 drops of this solution to about 2 cm<sup>3</sup> of solution E in a test-tube.

OBSERVATIONS

INFERENCES

(1 mark)

(1 mark)

- (iv) To about 2 cm<sup>3</sup> of solution E in a test tube, add 2 drops of aqueous potassium iodide.

OBSERVATIONS

INFERENCES

(1 mark)

(1 mark)

- (b) (i) Using a **metallic** spatula, ignite about one half of solid G in a Bunsen burner flame.

OBSERVATIONS

INFERENCES

(1 mark)

(1 mark)

- (ii) Place the other half of solid G into a boiling tube. Add 15 cm<sup>3</sup> of distilled water and shake well. Label this solution G. Use this solution for the following tests.

- I Place 2 cm<sup>3</sup> of solution G in a test-tube and determine its pH.

**OBSERVATIONS****INFERENCES**

(1 mark)

(1 mark)

- II To about 2 cm<sup>3</sup> of the solution obtained in (ii) above, add 3 drops of acidified potassium manganate (VII).

**OBSERVATIONS****INFERENCES**

(1 mark)

(1 mark)

- III To about 2 cm<sup>3</sup> of the solution obtained in (ii) above, add 2 drops of bromine water.

**OBSERVATIONS****INFERENCES**

(1 mark)

(1 mark)

- (iii) To the remaining solution G in the boiling tube, add the other half of solid F.

**OBSERVATIONS****INFERENCES**

(1 mark)

(1 mark)